

Prestan® AED Trainers

Prestan® Professional AED Trainer PLUS

The **Prestan Professional AED Trainer PLUS** provides quality and durability like no other trainer in the market. Its size and weight is comparable to an actual AED for a realistic training experience. The AED Trainer is compatible with a remote control.



Patent Pending
Unit Dimensions:
9.5" x 7.25" x 2.75"
Weight is 2.2 lbs.

PRESTAN PROFESSIONAL AED TRAINER PLUS

- PP-AEDT2-101 (Single AED Trainer, English/Spanish)
- PP-AEDT2-102 (Single AED Trainer, English/French)
- PP-AEDT2-401 (4-Pack AED Trainers, English/Spanish)
- PP-AEDT2-402 (4-Pack AED Trainers, English/French)
- PP-APAD2-1 (Single Set of Adult/Child Trainer Pads)
- PP-APAD2-4 (4-Pack Sets of Adult/Child Trainer Pads)

Additional languages available based on distributor-represented geographic area.

Optional Remote Control

PP-AEDT2-100-R

- Simulate loose pad
- Pause, play or stop
- Volume, language and scenario control
- Initiate audible "Press Deeper" during CPR
- Able to control multiple units at once



Prestan® AED UltraTrainer™

The **Prestan AED UltraTrainer** is a compact and customizable unit providing quality, durability and affordability. Sure to be an industry favorite.

PRESTAN AED ULTRATRAINER

- PP-AEDUT-101 (Single AED Trainer, English/Spanish)
- PP-AEDUT-102 (Single AED Trainer, English/French)
- PP-AEDUT-401 (4-Pack AED Trainers, English/Spanish)
- PP-AEDUT-402 (4-Pack AED Trainers, English/French)
- PP-UTPAD-1 (Single Set of dual Adult/Child Pads)
- PP-UTPAD-4 (4-Pack Sets of dual Adult/Child Pads)

Additional languages available based on distributor-represented geographic area.

Patent Pending
Unit Dimensions:
7.25" x 5.25" x 2.0"
Weight is 1.2 lbs.





Prestan® AED TrainerS - features

PRESTAN AED TRAINERS CAN BE CUSTOMIZED FOR ANY TRAINING ENVIRONMENT

The Professional AED Trainer PLUS and the AED UltraTrainer provide the options instructors have been waiting for, combining the most recent technology with realistic features not found on other AED Trainers.

Both AED Trainers are Equipped with Outstanding Features

- Clear, confident voice instructions that provide easy-to-follow directions
- CPR feedback/coaching
- Child defibrillation option
- Training pads are pre-connected and ready to use, unlike other AED trainers that use outdated plug-in technology
- Pre-configured shock scenarios (5 options available)
- Semi-or-fully automatic shocking option
- Compression metronome option (turn on / turn off)
- Rescue Breaths option (turn on / turn off)
- Unit programming can easily toggle between 2 languages
- Programmed with the most up-to-date software and language capability in full compliance with current industry guidelines
- Updates are simple with a module (Professional AED Trainer PLUS) or dongle (AED UltraTrainer)
- Prestan AED Trainers carry a three (3) year warranty

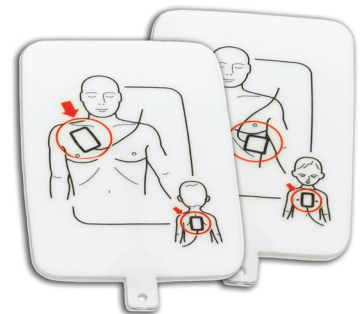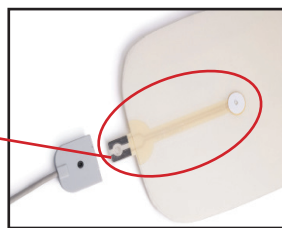

Both AED Trainers are Equipped with Amazing Training Pads

- Pads have an embedded 'switch' which senses proper pad attachment to a manikin. This **Pad Sensing System** (patent pending) allows the Prestan AED Trainers to be used with almost any CPR manikin - no external device needed
- The silicone base adhesive on the pads, coupled with the foam structure, ensures easy application & removal with no residue
- Pads can be used approximately 100 times (with proper care) during CPR/AED classes
- Dual-graphics pads work for both adult and child training